


Connection Instruction

Step 1, Checking the operating system of computer, make sure it is 32-bit or 64-bit, install corresponding driver software

Step 2, Turn on instrument, connect it to the computer through data line, right-click “my Computer”→Management→Device Manager→Ports COM and LPT, you can view the corresponding port

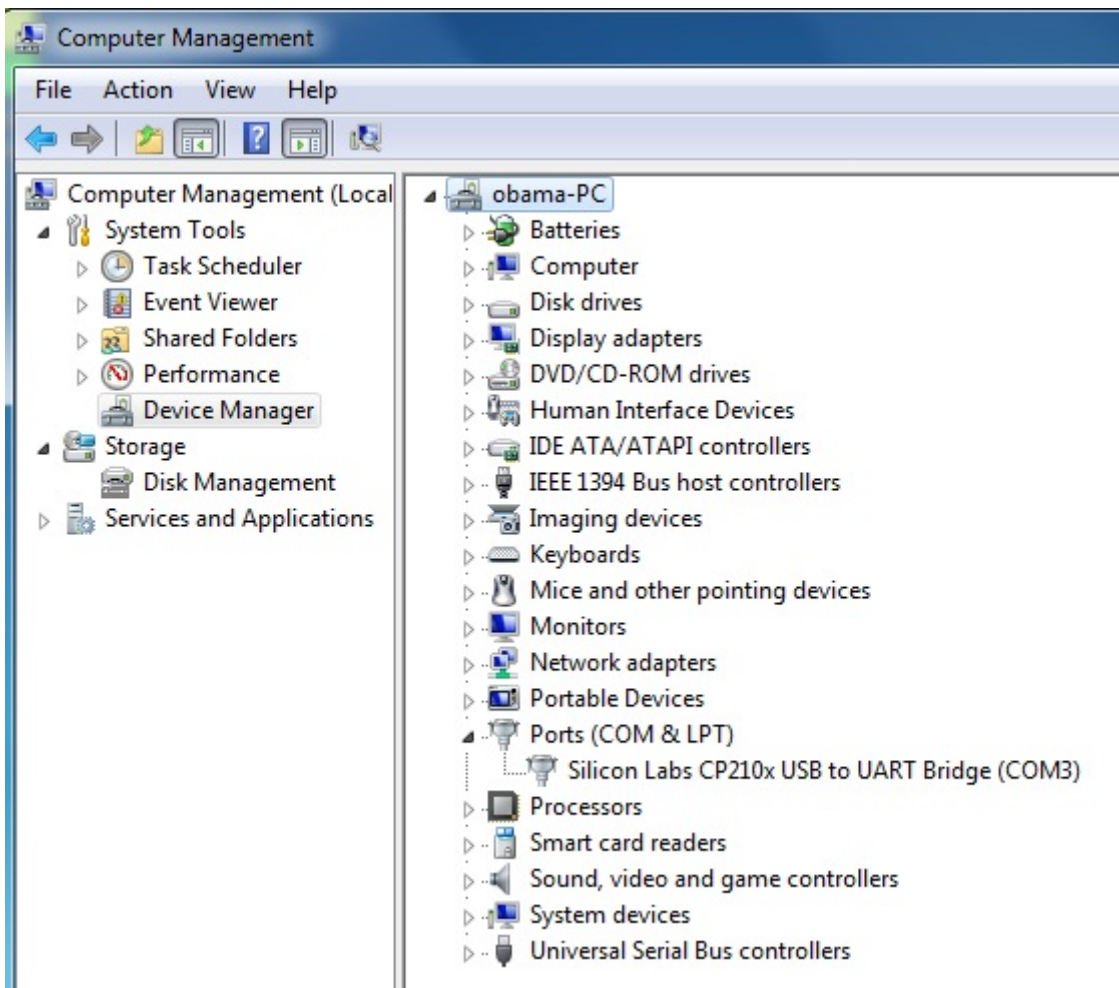


Note: If the port number exceeds COM30, you need to modify it. Here is the steps

Right-click “my Computer”→Management→Device Manager→Ports COM and LPT, you can view the corresponding port, choose COM1→Right-click→Property→Port Settings→Advanced→Choose the right port among COM1~COM30 in Pull-down Menu of Window “COM port number” → Confirm, exit to “Device Manager” → Choose the original port → Right-click→ Scan for hardware changes. Then port number has been changed.
Operation interface as below:



Step 3, Check out the corresponding port, open the digital software , select the corresponding port, click “open”→click “read” to import data, then you can view the saved data. If click “excel”, the data will be saved into excel spreadsheet and then can be printed. If click “delete”, the data will be deleted. (Recommend to save the data into excel document, and then delete the data in the instrument daily.)

Operation interface as below:

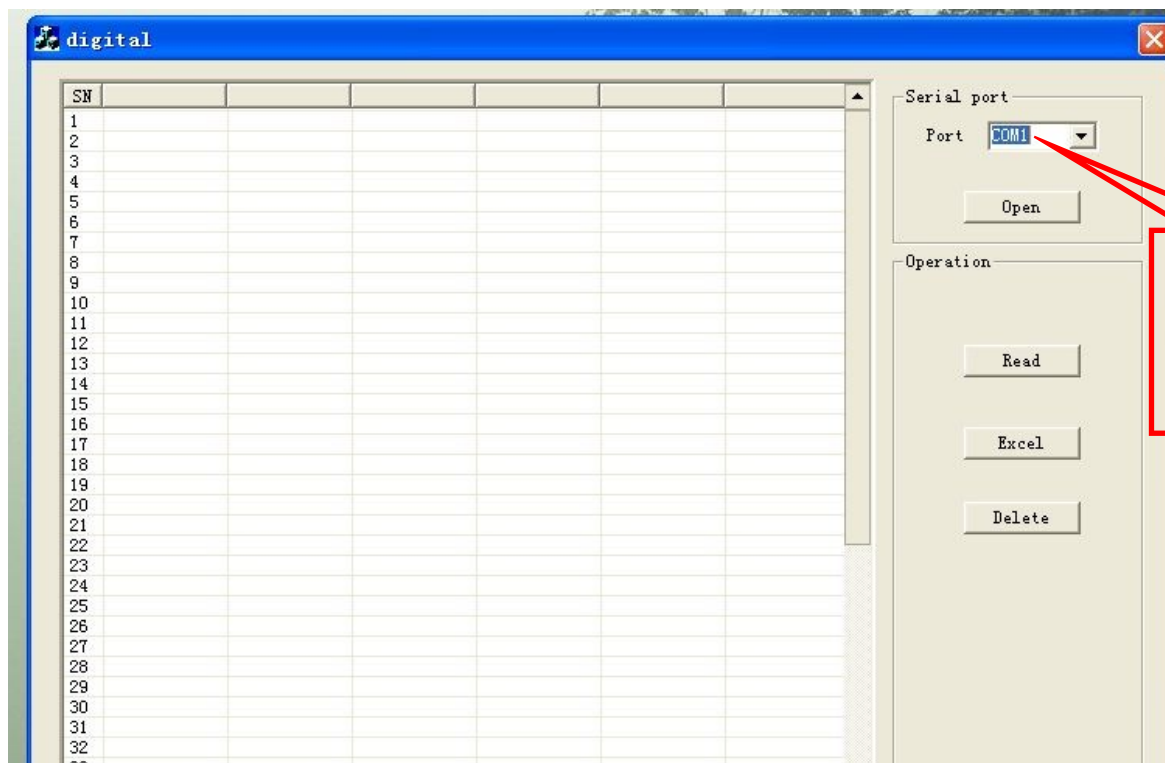

Choose the port referred in Step 2

